

FastRuler Middle Range DNA Ladder, ready-to-use

Catalog Number SM1113

Pub. No. MAN0013028 Rev. D00



WARNING! Read the Safety Data Sheets (SDSs) and follow the handling instructions. Wear appropriate protective eyewear, clothing, and gloves. Safety Data Sheets (SDSs) are available from [thermofisher.com/support](https://www.thermofisher.com/support).

Contents and storage

Cat. No.	Contents	Amount	Storage
SM1113	FastRuler Middle Range DNA Ladder, ready-to-use	2 x 500 µL	at room temperature or at 4 °C for periods up to 24 months.*
	6X MassRuler DNA Loading Dye	1 mL	

*For longer periods store at -20 °C.

Description

Thermo Scientific™ FastRuler™ Middle Range DNA Ladder is specially designed for fast sizing and quantification of double-stranded DNA in 48-well (or 96-well) gels as well as in conventional agarose gels.

The Ladder consists of five blunt-end chromatography purified individual DNA fragments (in base pairs): 5000, 2000, 850, 400 and 100, mixed in equal quantities. The ladder bands are easily resolved after an 8-14 min in run on appropriate agarose gels in a short 10-20 mm linear separation distance and visualized by ethidium bromide or SYBR™ Green I staining.

The Ladder is premixed with a loading buffer and can be applied directly onto an agarose gel.

Storage and Loading Buffer

10 mM Tris-HCl (pH 7.6), 10 mM EDTA, 0.005 % bromophenol blue and 10 % glycerol.

6X Thermo Scientific™ MassRuler™ DNA Loading Dye

10 mM Tris-HCl (pH 7.6), 0.03 % bromophenol blue, 60 % glycerol and 60 mM EDTA.

Protocol for Loading


Step 1: Mix gently

Step 2: Load 1 µL per 1 mm gel lane

Recommendations

- Do not heat before loading.
- For sizing:
 - dilute your DNA sample with the MassRuler DNA Loading Dye (#R0621, supplied with the ladder): mix 1 volume of the dye with 5 volumes of the DNA sample;
 - load equal volumes of the Ladder and an experimental sample;
- For quantification:
 - adjust the concentration of your sample such that the expected amount of material loaded is approximately equal to that of a DNA Ladder band of a nearest size.
- For DNA band visualization with SYBR™ Green and other intercalating dyes, do not add the dyes into the sample, use gel staining after electrophoresis or include dyes into agarose gel to avoid aberrant DNA migration.
- **Important note:** For DNA bands visualization with GelRed™ use gel staining after electrophoresis to avoid aberrant DNA migration.

Amount of DNA (ng) in each band of FastRuler Middle Range DNA Ladder, ready-to-use

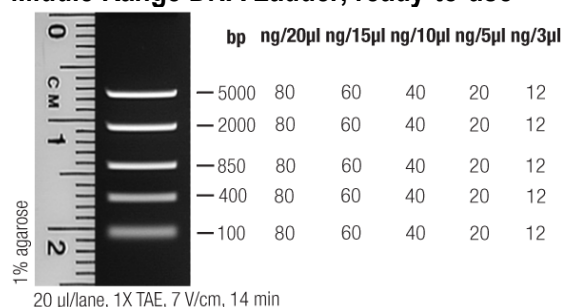

Note

The apparent intensity of bands containing equal ng quantities of DNA may differ in different horizontal sections of gel (diminishes from top to bottom). This effect derives from higher diffusion rates of the shorter fragments, resulting in the lower compactness and lower peak intensity of the corresponding bands. Among factors affecting the severity of this effect are: temperature (the lower the better), overall duration of an electrophoretic and staining procedure (shorter), percentage of gel (higher). The time course of DNA fragments separation is presented in the Fig.1.

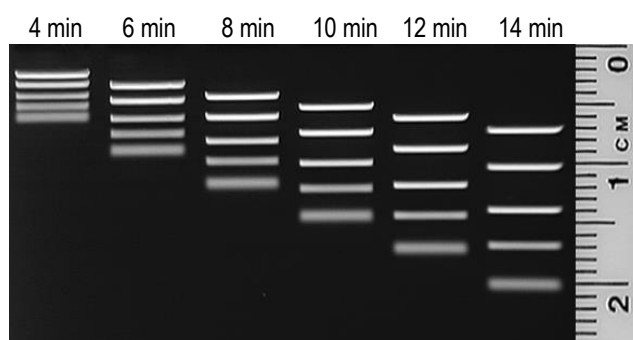


Fig. 1. Time course of band separation.
Electrophoresis conditions: 1 % agarose, 1X TAE, 7 V/cm.

Revision history: Pub. No. MAN0013028

Revision	Date	Description
D00	2025-04-08	Extended storage at room temperature to 24 months

Limited product warranty

Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale at www.thermofisher.com/us/en/home/global/terms-and-conditions.html. If you have any questions, please contact Life Technologies at www.thermofisher.com/support.



Thermo Fisher Scientific Baltics UAB | V.A. Graiciuno 8, LT-02241 Vilnius, Lithuania
For descriptions of symbols on product labels or product documents, go to thermofisher.com/symbols-definition.

The information in this guide is subject to change without notice.

DISCLAIMER: TO THE EXTENT ALLOWED BY LAW, THERMO FISHER SCIENTIFIC INC. AND/OR ITS AFFILIATE(S) WILL NOT BE LIABLE FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE, OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING YOUR USE OF IT.

Important Licensing Information: These products may be covered by one or more Limited Use Label Licenses. By use of this product, you accept the terms and conditions of all applicable Limited Use Label Licenses.

©2025 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified. GelRed is a registered trademark of Biotium Inc.

thermofisher.com/support | thermofisher.com/askaquestion

thermofisher.com

ThermoFisher
SCIENTIFIC